



Series	BT300 Series 300 µl Barrier Tip	2090 Series 300 µl Universal Tip		
Part Number	BT300	2090 2090.N	2097 2097.N 2097.NS 2097.S	2092 2092.N 2092.NS 2092.S
Graduation Marks	Graduation marks are indicated at both 0.5 µl, 2 µl and 10 µl locations Refer to product image for visuals			
Tip Composition	Neptune pipette tips are made of virgin polypropylene			
Tip Type	S ³ - Low Retention			
Filter Material	High Density Polyethylene Filter	Non-Filtered Products		
Offered in Sterile Format	Yes	No	Yes	
Configuration	Racked	Bulk	ESP Reload	Racked
Packaging Breakdown	96 tips/ rack 10 racks/ pack 5 packs/ case	1000 tips/ bag 20 bags/ case	96 tips/ card 20 cards/ pack 4 packs/ case	96 tips/ card 10 racks/ pack 5 packs/ case
Autoclavable	No	Autoclavable at 120°C for 10-15 minutes at 15 PSI		
Storage Condition	Store in a clean, dry environment at room temperature 15-30°C			

BT300 Series

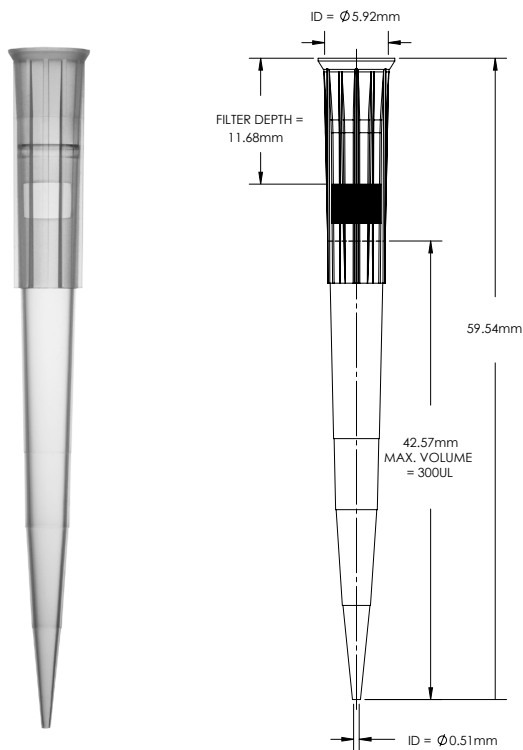

2090 Series

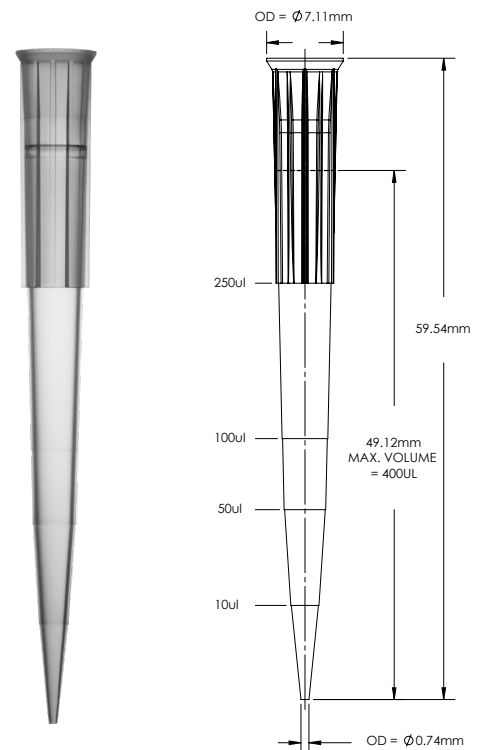

Product Specification Sheet

Quality Control:

Certificates of Compliance	Each lot undergoes stringent inspection and individual lot testing ensures Neptune products are certified RNase, DNase, DNA and Endotoxin-free. Visit www.neptunescientific.com to obtain a copy of a certificate of compliance for your Neptune product.
RNase/ DNase	Products are washed in distilled water and concentrated via centrifugation. Samples are added to previously established nucleic acid standards, incubated for one hour at 37°C, and tested on a 2% gel using electrophoresis. Products must show no degradation of standards to pass. Test sensitivity is 10-7 Kunitz units/µl.
Nucleic Acid	Products are washed in distilled water and concentrated via centrifugation. Then, samples are added to protocol specified PCR reactions and thermal cycled for 50 cycles. A 2% agarose gel electrophoresis is used to examine experimental samples, positive controls, and negative controls. To pass, product samples must show no DNA amplification. Test sensitivity is 10 ng.
Endotoxin/ Pyrogen	Products are tested for endotoxins by using the Limulus Amebocyte Lysate (LAL) gel assay according to FDA guidelines. Test sensitivity is 0.06 EU/ml.
Sterilization	Products are sterilized using electron beam irradiation.
Traceability	Each product contains a 5 digit lot number located on the rack, pack and case of each finished good. With Neptune's advanced manufacturing process all raw materials are able to be traced for maximum quality assurance.

Advancements in Liquid Handling

S³	Neptune's exclusive S ³ polymer was designed to increase pipetting accuracy by virtually eliminating tip retention and sample hold-up.
ESP Reload	Neptune's ESP (Environmental Sustainable Pack) was the industry's first pipette reload system designed to minimize plastic waste by 90% and provide an environmentally friendly solution.
Aerosol Barrier Tip	Specifically engineered to reduce cross contamination.

Pipettor Compatibility:

Biohit™ M100 and M200

Biohit Proline Plus™ 100 µl
(works with 2090 Series)

Biohit Proline Plus™ 200 µl

Brand Transferpette S™ 200 µl

Brand Transferpette Electronic™ 300 µl

Capp™ 50 µl, 100 µl and 300 µl

CLP Beta-Pette™ 200 µl

CLP Poseidon™ 50 µl, 100 µl, 200 µl
and 300 µl

CLP Poseidon Electronic™ 200 µl

Eppendorf Reference™ 200 µl

Eppendorf Research™ 100 µl
(works with 2090 Series)

Eppendorf Research™ 200 µl and 300 µl

Eppendorf Research Plus™ 200 µl

Eppendorf Xplorer™ 300 µl

Finnpipette™ 50 µl, 200 µl and 300 µl

Finnpipette™ Electronic 300 µl

Gilson Pipetman™ P200

Gilson Pipetman Ultra™ U200

Hamilton™ 25 µl, 100 µl and 300 µl

Nichiryo Nichipet EX™ 200 µl

Nichiryo Oxford Multimate™ 300 µl

Socorex Calibri 822™ 200 µl

VWR Ergonomic High Performance™ 20 µl
(works with 2090 Series)

VWR Ergonomic High Performance™ 200 µl

VWR Ultra High Performance™ 20 µl,
100 µl (works with 2090 Series)

VWR Ultra High Performance™ 200 µl